


ACCU-Beam®

TTI Medical

ACCU-Beam® Universal CO₂ Laser Laparoscope Coupler....

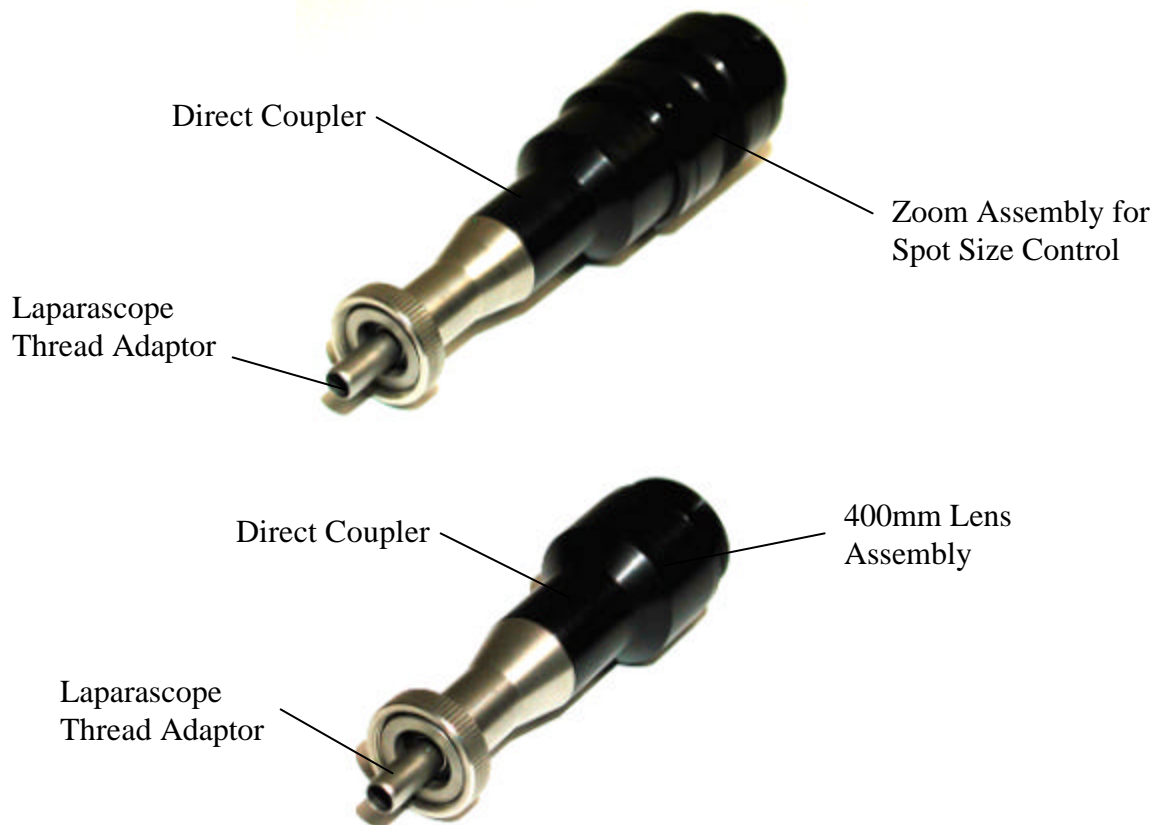

The ACCU-Beam® Universal CO₂ Laser Laparoscope Coupler directly attaches to the proximal end of the operating laparoscope channel. The alignment stem aids in beam centration through the operating channel. Laparoscope Thread Adaptors with alignment stems are available for all commonly used laparoscopes.

The coupler can be fitted with either a 400mm focusing lens or a 250mm to 450mm zoom focusing assembly. The zoom focus assembly enables beam focusing and defocusing while keeping the laparoscope stationary.

The 400mm focusing lens and zoom focusing assembly fits directly to articulated laser arms with 15/16-12 threads. There is a wide selection of thread adaptors for arms with different thread specifications.

ACCU-Beam® Universal CO₂ Laparoscope Coupler

CO₂ Laser Laparoscope Coupler

3113	Direct Coupler
3116	Zoom Focusing System 180-400mm
3125	400mm Lens Assembly

Laparoscope Thread Adaptors

3181AS6	Standard - for Storz, Wisap, Figert, Leisegang and Gynoscope - 6mm
3182AS5	Olympus - 5mm
3183AS6	Cabot - 6mm
3184AS5	Wolf - 5mm
3184AS8	Wolf - 8mm
3185AS7.5	Storz Bayonet - 7.5mm

CO₂ Laser Articulated Arm Thread Adaptors

1101	Sharplan Quick Disconnect - for models 1040, 1060 and 1100
1102	LaserSonics / NIIC
1103	LaserSonics / Merrimack Lumina 40 Silver Arm / Paragon 50
1104	Coherent / Xanar
1105	Coherent XA-30
1106	Coherent 450 & 451
1108	Zeiss
1109	Sharplan Twist-Lock - for models 1020, 1050, 1055, 1075 and Ultra Pulse
1110	Union Medical
1111	Lasering
1112	LaserSonics LS-500
1115	Coherent XA-50

Note: Sharplan 1060 has two arm versions - 1101 & 1109. Thread adaptors are not required for Surgilase, Laser Engineering, Sharplan 720, 733A, 734, and LaserSonics Illumina 40 Black Arm, and Dornier 401.

ISO-9001:1994
AMTAC Cert.No. 163



TRANSAMERICAN TECHNOLOGIES INTERNATIONAL
2246 Camino Ramon • San Ramon, CA 94583 • USA
Tel: 925-355-0750 • 800-322-7373 • Fax: 925-355-0777
e-mail: info@ttimedical.com • website: www.ttimedical.com

Rev060903