


ACCU-Beam®

TTI Medical

ACCU-Beam® Universal Zoom Video EndoAdaptor for Rigid and Flexible Endoscopes....

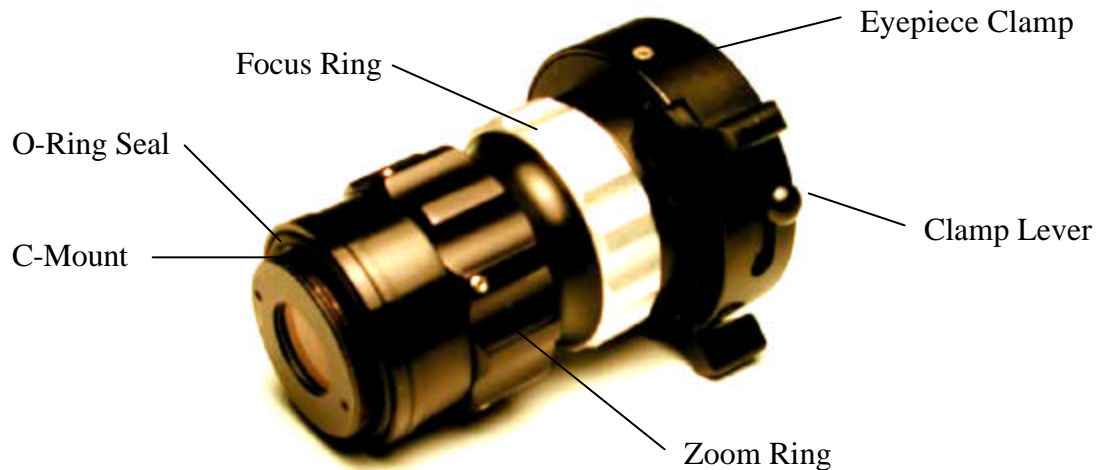

U.S. Patent No. 6,113,533 and 6,115,973

The ACCU-Beam® Universal Zoom EndoAdaptor features a proprietary parfocal optical design. Image magnification can be continuously changed through a focal range of 13mm to 28mm without the need to adjust focus. Image quality remains sharp and clear.

A C-mount thread is used for attaching standard C-mount video cameras. The Zoom EndoAdaptor can be modified to accept other camera types.

The eyepiece clamp accommodates a standard 32mm eyepiece. The eyepiece can be locked in place or unlocked to rotate the scope within the clamp. Eyepiece adaptors are available for a wide variety of flexible endoscopes.

ACCU-Beam® Universal Zoom EndoAdaptor

P/N 8019

Specifications:

Zoom Range	Parfocal from 13mm to 28mm
Focus	One time adjustment – video image remains in focus through entire zoom range.
Camera Mount	C-mount and custom camera mounts.
Eyepiece Clamp	For 32mm rigid scope eyepieces – lockable.
Soak Disinfecting	Use 2% activated glutaraldehyde.
Sterilization ETO, Steris or Sterrad	
Length	2.5 inches
Weight	4.3 oz

Flexible Endoscope Adaptors: Connects flexible scope eyepiece to eyepiece clamp.



8050 Pentax	8051 Machida
8052 Olympus OES	8053 Olympus pre-OES
8054 Fujinon	8055 Schott LS7
8056 Schott FS100 & FS101	

ISO9001:1994
AMTAC Cert.No.163



TRANSAMERICAN TECHNOLOGIES INTERNATIONAL
2246 Camino Ramon • San Ramon, CA 94583 • USA
Tel: 925-355-0750 • 800-322-7373 • Fax: 925-355-0777
e-mail: info@ttimedical.com • website: www.ttimedical.com

Rev060903